


Thank you ECL™ Plus, you've been great.

For every success there's a successor. It's called progress. So while we knew we had something special in Amersham™ ECL Plus, our Western blotting team was quietly working on the next generation of detection reagent. The result: a new substrate that operates with superior levels of sensitivity, signal intensity and stability than even its famous predecessor, making it an excellent choice for CCD imagers. **Welcome to ECL Prime.**



Find out more about Amersham ECL Prime at
www.gelifesciences.com/eclprime

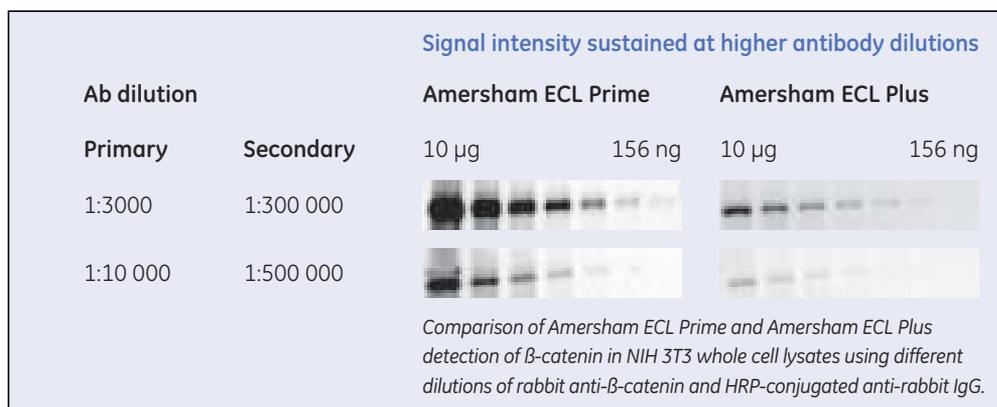

imagination at work



Learn how current events are impacting your work.

ScienceInsider, the new policy blog from the journal *Science*, is your source for breaking news and instant analysis from the nexus of politics and science.

Produced by an international team of science journalists, *ScienceInsider* offers hard-hitting coverage on a range of issues including climate change, bioterrorism, research funding, and more.

Before research happens at the bench, science policy is formulated in the halls of government. Make sure you understand how current events are impacting your work. Read *ScienceInsider* today.

www.ScienceInsider.org

ScienceInsider

Breaking news and analysis from the world of science policy



MPC-200 Multi-manipulator system

Versatile: User friendly interface controls up to two manipulators with one controller. Select components to tailor a system to fit your needs.

Expandable: Daisy chain a second controller and operate up to four manipulators with one input device.

Stable: Stepper motors and cross-rolled bearings guarantee reliable, drift-free stability.

Doubly Quiet: Linear stepper-motor drive reduces electrical noise. Thermostatically-controlled cooling fans barely whisper.

Make the right move!



SUTTER INSTRUMENT

PHONE: 415.883.0128 | FAX: 415.883.0572
EMAIL: INFO@SUTTER.COM | WWW.SUTTER.COM

1,000s OF GRANTS MILLIONS IN FUNDING



GrantsNet. The first comprehensive science grants database.

GrantsNet is expanding its listings of some 900 funding programs from private foundations and not-for-profit organizations to include 400 to 500 new entries from the grants.gov site. **This provides the first comprehensive database of funding opportunities** to research scientists and administrators, career counselors, financial aid specialists, and undergraduate and graduate students. For listings, go to www.grantsnet.org

Science Careers

From the journal *Science*

INTRODUCING

AAAS MemberCentral



The exclusive new website for the AAAS member community.

AAAS MemberCentral is a new website focused on helping you — the scientists, engineers, educators, students, policymakers, and concerned citizens who make up the AAAS community — connect like never before.

On MemberCentral you can contribute to discussion groups or blogs, participate in a webinar, or share photos of your field research. You can exchange ideas, learn about your fellow members, and gain fresh insights into issues that matter to you the most. MemberCentral is also an easy access point for a wide variety of other AAAS membership benefits, like discounts on cars and books, travel opportunities, and more.

Experience MemberCentral for yourself. Visit MemberCentral.aaas.org today and log in using your *Science* online username and password.



MemberCentral.aaas.org

HPLC COLUMNS

The new Synchronis range of high-performance liquid chromatography (HPLC) columns unerringly delivers consistent, predictable separations, from run to run and column to column. When developing a new method, one of the key goals for the chromatographer is to achieve a consistent, reproducible separation. The selection of a highly reproducible HPLC column is essential to attaining this goal. The new Synchronis column range has been engineered to deliver exceptional reproducibility by providing highly pure, high surface-area silica, dense bonding, and double endcapping, all controlled and characterized through rigorous testing. New enhanced, automated column packing methods used within the Synchronis range deliver greater consistency, and every column is individually tested to ensure that it meets the required specifications. These extensive testing and quality control procedures ensure the delivery of a completely consistent product, column after column. Synchronis columns are available with 5 μm particle size for conventional HPLC applications and with 1.7 μm particle size for high-speed, high-efficiency UHPLC separations.

Thermo Fischer Scientific

For info: 800-532-4752 | www.thermoscientific.com/Synchronis



PRESSURE/FLOW RATE SENSORS

The new Mitos Sensor Units provide a flexible system for measuring and displaying pressure and flow rates in microfluidic systems. With a real-time display and low internal volumes, each sensor is designed to minimize interference with the liquid flow. All data is clearly displayed on the Mitos Sensor Display unit, which interfaces with each of the sensor units using a simple push-and-click action. Comprising five flow rate sensors and one pressure sensor (-0.5–30 bar), the complete range is fully interchangeable and each device can run stand-alone or integrated with a Mitos P-Pump for the efficient logging of flow rate output. All data is easily transferrable to a PC for further analysis via USB. Ongoing enhancements to the system will soon enable the flow rate sensors to communicate with the Mitos P-Pump to create a closed loop pumping system. This will enable the user to automatically control flow rates using pressure, providing the ultimate in smooth liquid flow.

Dolomite

For info: +44-1763-242491 | www.dolomite-microfluidics.com

CONTROLLED RATE FREEZER

The Bio-Cool is the only controlled rate freezer that does not require expendable liquid nitrogen and the associated pumping, refilling, and storage challenges posed by a cryogenic liquid. The Bio-Cool, simply plugs into a standard electrical outlet and quietly provides low temperature cooling to either -40°C or -80°C. Unlike liquid nitrogen systems, which surround samples with a cold vapor phase, the Bio-Cool immerses samples into a well-circulated, cold liquid bath. The circulating liquid allows more efficient heat transfer to the samples and maintains a more consistent temperature profile at all locations. The Bio-Cool comes standard with an easy-to-use programmable microprocessor controller that enables programming

and storage of up to eight different protocols, each with up to eight distinct segments specifying ramp rate, hold temperature, and hold time. Elimination of liquid nitrogen handling and usage, better freezing results, simple straightforward operation, and world-class service support all make the Bio-Cool the controlled rate freezer of choice for most clinical, commercial, and research applications.

SP Scientific

For info: 845-255-5000 | www.spscientific.com

IMAGING CAMERA

Engineered for high-sensitivity and high-speed performance in low-light imaging applications, the Rolera EM-C2 camera aids researchers in applications such as spinning-disc confocal imaging, total internal reflection fluorescence (TIRF) microscopy, ratiometric ion imaging, and fluorescence recovery after photobleaching (FRAP). Fast frame rates essential to cutting-edge live cell imaging studies are enabled by the camera's 40 MHz pixel clock rate. Up to 34.2 full-resolution (one megapixel) frames per second can be read out and transferred over an optimized 800 Mb/s implementation of the IEEE 1394b FireWire protocol. The Rolera EM-C2 delivers very low noise and high-sensitivity across a broad spectrum common to fluorescence experiments. This camera introduces an Easy-EM mode, which optimizes camera EM Gain settings with a single click. Researchers no longer have to gauge sensitivity using the arbitrary sliders that are required with the use of existing cameras to estimate necessary EM Gain. The camera includes the QCapture Suite software, which can be installed and operational within minutes.

QImaging

For info: 800-874-9789 | www.qimaging.com

Electronically submit your new product description or product literature information! Go to www.sciencemag.org/products/newproducts.dtl for more information. Newly offered instrumentation, apparatus, and laboratory materials of interest to researchers in all disciplines in academic, industrial, and governmental organizations are featured in this space. Emphasis is given to purpose, chief characteristics, and availability of products and materials. Endorsement by *Science* or AAAS of any products or materials mentioned is not implied. Additional information may be obtained from the manufacturer or supplier.

Science Without Borders

17–21 February, Washington, D.C.

AAAS, publisher of *Science*,
thanks the sponsors and supporters of
the *2011 Annual Meeting*

Presenting sponsor



SUBARU

Canada



EUROPEAN COMMISSION
European Research Area



THE  KAVLI FOUNDATION



To become an exhibitor or
sponsor contact:

Jill C. Perla
Manager of Marketing, Exhibits,
and Sponsors
E-mail: jperla@aaas.org
Phone: 202-326-6736

Web site:
<http://www.aaasmeeting.org>

“How do we know
this lead molecule
is novel?”

“SciFinder—
of course.”

Need to assess the novelty of substances?

SciFinder is the answer.

It includes CAS REGISTRYSM the most comprehensive substance information available, integrated with relevant journal articles and patents.

Give your research team the highest quality and most timely scientific information resource.

Make SciFinder an essential part of your research process.

For more information about SciFinder, visit www.cas.org or e-mail help@cas.org.

an essential
✓
SciFinder®—Part of the process.™



CAS is a division of the American Chemical Society

www.cas.org